

SURGICAL OFFICE SET-UP: INSTRUMENTS AND EQUIPMENT



KEN GROSS , M.D.

CONFLICTS OF INTEREST for KEN GROSS M.D.

- 3M Pharmaceuticals/Graceway
- Abbott Pharmaceuticals
- Activis Pharmaceuticals
- Allergan Inc.
- Altana Inc.
- Atrix Laboratories, Inc.
- Azur Tech, Inc.
- Barrier Pharmaceuticals
- Bertek Pharmaceuticals Inc.
- Bristol-Myers Squibb
- Celgene Corporation
- Cipher Pharmaceuticals
- Collagenex/Galderma
- Connetics Corporation
- Connective Therapeutics
- Corixa Corporation
- Cytotech Pharmaceuticals
- Dermik Laboratories Inc.
- Derm Tech
- Dow Pharmaceuticals
- Dusa Pharmaceuticals
- Electro Optical Sciences
- Fujisawa/Astellas
- Genentech, Inc. // Roche
- Genomics
- Healthpoint Laboratories
- Hill Dermaceuticals
- Hyal Pharmaceuticals
- Leo Pharma
- Matrix Pharmaceuticals
- Meiji
- Novartis Pharmaceuticals
- Obagi
- Organogenesis Inc.
- Paddock Laboratories
- Parke-Davis Pharmaceuticals
- Leo Pharma (Picato)
- Pennederma Inc.
- Pfizer
- Procyon Pharmaceuticals
- Promius Pharma
- Quintiles Laboratories
- Sandoz Pharmaceuticals
- Teva Pharmaceuticals
- Tolmar, Inc.
- Valeant Pharmaceuticals

I have done Clinical Studies or Lectured for these companies

To add your name or company, call Skin Surgery Medical Group Inc. at (858) 292-5101

Additional **CONFLICT OF INTEREST:**

I would like to thank **Brian Kelley** from **George Tiemann & Co.** for helping me obtain photos of surgical instruments used in this lecture and providing instruments for the closure lab.

He has been a great help to both the “Fundamentals of Mohs Surgery” and the “Closure” Courses.

I receive **no special discounts or consideration from George Tiemann & Co.**, other than their usual great service.

The nature of this talk requires that I recommend specific products by name or manufacturer.

I have no financial interests in any product or company mentioned in this lecture.

The nature of this talk requires that I recommend specific products by name or manufacturer.

I have **no financial interests** in any product or company mentioned in this lecture.

Do Not Distribute or Duplicate



should not be a
to the grave with
the intention of
arriving safely in an
and well preserved
body, but rather
dead in swimmy,
in one hand,
Lotto
in the other,
body thoroughly
used up, totally worn
out and screaming
"UO HOO"
What a ride!

Payment Due at Time
of Cosmetic Services

Hand Sanitizer for Patients

Do Not Distribute or Duplicate



*should not be a
to the grave with
the intention of
arriving safely in an
and well preserved
body, but rather
to die in swimming,
in one hand,
Lotto
in the other,
body thoroughly
used up, totally worn
out and screaming
"UO HOO
What a ride"*

Operating Room Accreditation

Payment Due at Time
of Cosmetic Services

Attention Patients:
Effective April 14, 2011
We'll have a special package for the
before the procedure:
\$27.00
\$34.00
\$39.00
\$43.00
\$49.00
\$57.00
\$64.00
\$74.00
\$84.00
\$12.00

Hand Sanitizer for Patients



INSTITUTE FOR MEDICAL QUALITY
AMBULATORY PROGRAM
A subsidiary of the California Medical Association



Accreditation

SKIN SURGERY MEDICAL GROUP (10010)
5222 Balboa Avenue, Sixth Floor
San Diego, CA 92117

*has completed all requirements for accreditation as an ambulatory care facility,
and has demonstrated compliance with the standards
of the IMQ Ambulatory Program for the period*

July 29, 1999 - July 28, 2002

KM Tan, MD
Chair, IMQ Board of Directors

Robert M. Takamoto, MD
Chair, Ambulatory Care Review Committee

Do Not Distribute or Duplicate



Oversized Doorway on Magnetic "Keep Open"

Pocket Doors are Efficient:
Use ONLY Commercial Grade Double Hung



Do Not Distribute or Duplicate



Pocket Doors are Efficient:
Use Double Hung Only



“High Speed” Mirrors



Do Not Distribute or Duplicate




ACLS



Elman "Chair Ground" no longer available

Sterile Tip For Each Patient

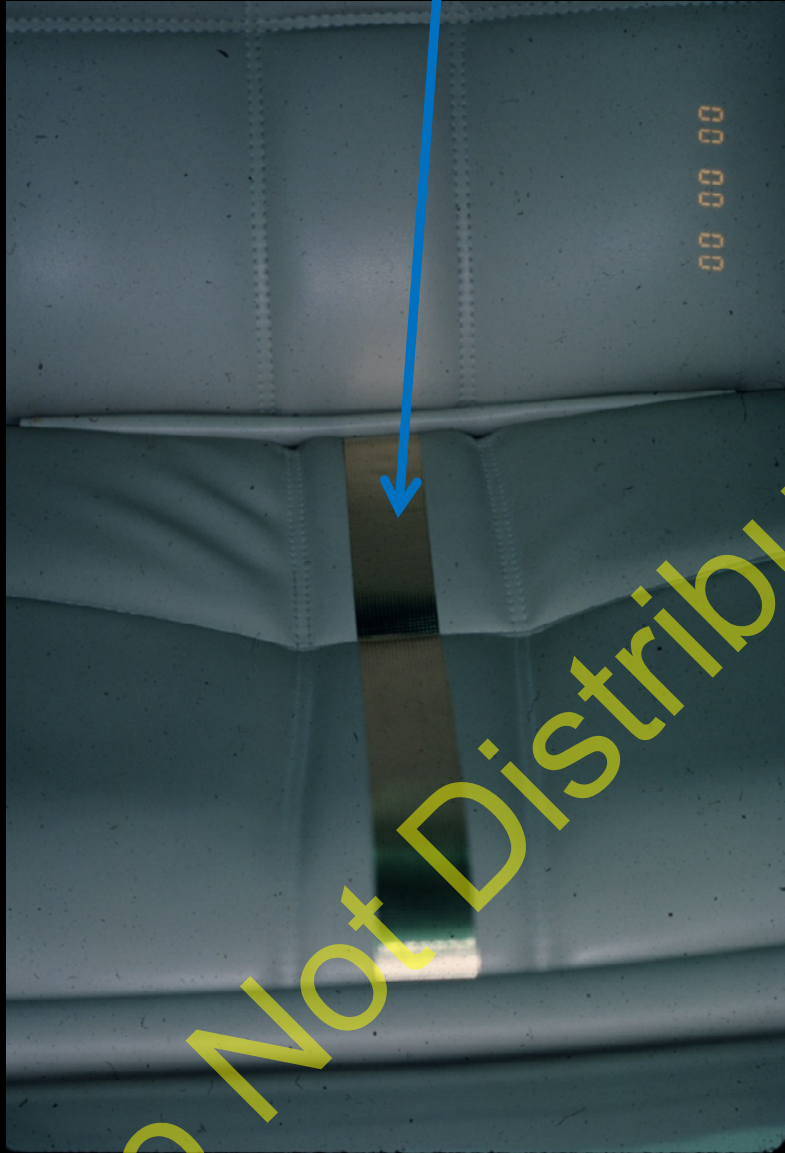

Do Not Distribute Or Duplicate



More Char and More Smoke and Less ability for Pin Point Hemostasis




Old Gomco

New Cheaper Suction



**Adjustable Chairs
on Both Sides of
the Table**

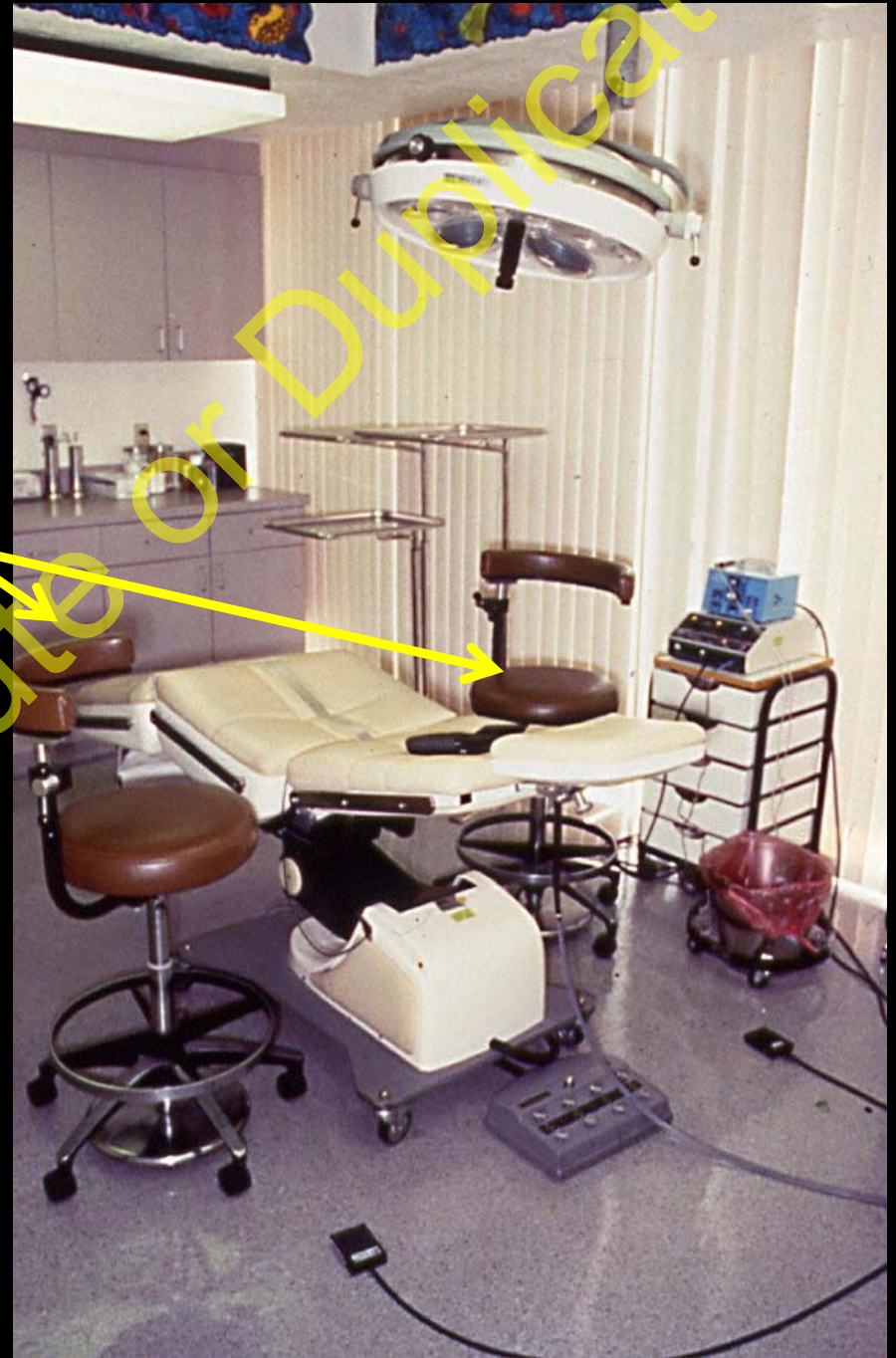

**Any Foot-Operated,
Adjustable O.R. Chair
Will Do**



Do Not Distribute or Duplicate

**Big Trays:
Sturdy,
Stackable**




Big Trays:
Sturdy,
Stackable

Easily
Adjustable



What Your Tray Should NEVER Look Like

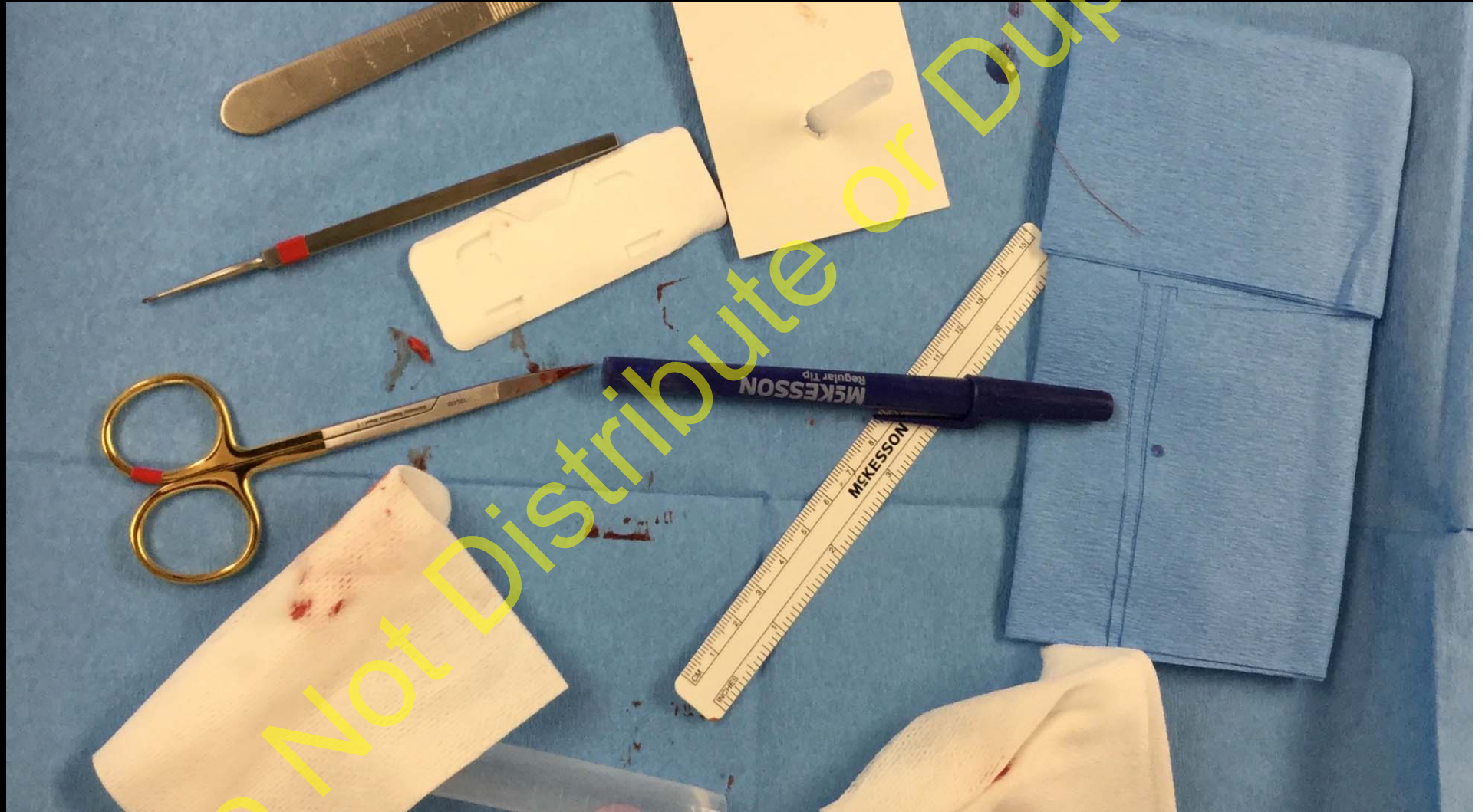


What Your Tray Should ALWAYS Look Like



Per OSHA, the only items that go in the red bags are items that are saturated with human bodily fluids like blood--- the best way to think of it is, if it is dripping or if it can be wrung out it's saturated and goes in the red bucket
RED BUCKET== \$\$\$

New OSHA Rule

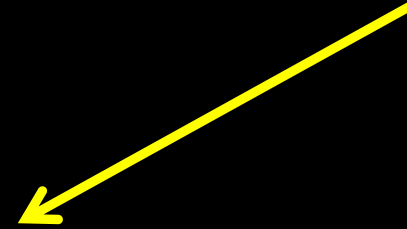

Do Not Distribute or Duplicate



Do Not Distribute or Duplicate

NEEDLES
STERILIZED
SHARP
Authorized
Personnel
Only
DANGER

Simply Saline
Wound Wash

ADHESIVE
Vi-Drape
by MCD
DANGER: EXTREMELY FLAMMABLE
WARNING: CONTENTS UNDER PRESSURE
(SEE OTHER CAUTIONS ON LABEL)
NET WT. 5oz (142g)
LOT NUMBER AND EXPIRATION DATE
Medical Device
MADE IN USA

Do Not Distribute or Duplicate

Vi-drape (Adapt) skin adhesive



Do Not Distribute or Duplicate



THINNER



THICKER


These Are The Most Comfortable Gloves I Have Used

I HAVE NO FINANCIAL INTEREST IN THIS COMPANY

Do Not Distribute or Duplicate

Surgical Gloves

- Surgical gloves

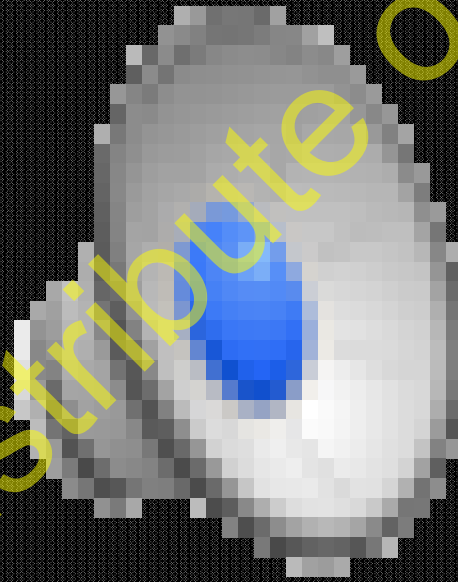

Do Not Distribute or Duplicate



Wet 4x4 In Ziploc Bag



Do Not Distribute or Duplicate



Wet 4x4 In Ziploc Bag

**We Make Ice Packs For
Patient Use in Office**



Do Not Distribute or Duplicate



Anesthetic is injected before full prep and gloving. This allows EPI to have time to work.



Do Not Distribute or Duplicate



Anesthetic injected before full prep and gloving up so EPI has time to work.

When sterile, we have extra single sterilized gloves to put over surgery gloves if more anesthesia is needed.



Do Not Distribute or Duplicate

Quick, Easy to Use, and "Hands-Off"



Blade Remover

I HAVE NO FINANCIAL
INTEREST IN THIS COMPANY

Quick, Easy to Use, and "Hands-Off"







Blade Remover

Sharps Disposal

I HAVE NO FINANCIAL
INTEREST IN THIS COMPANY

Monsel's
Exp:10/12

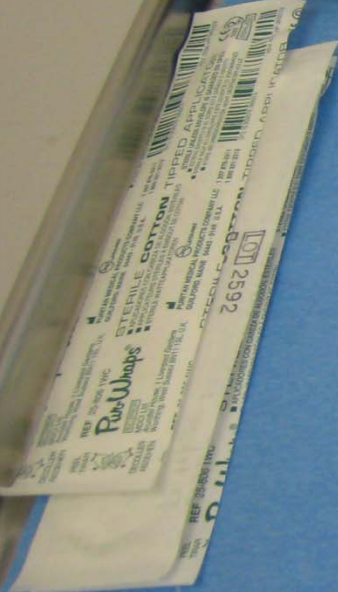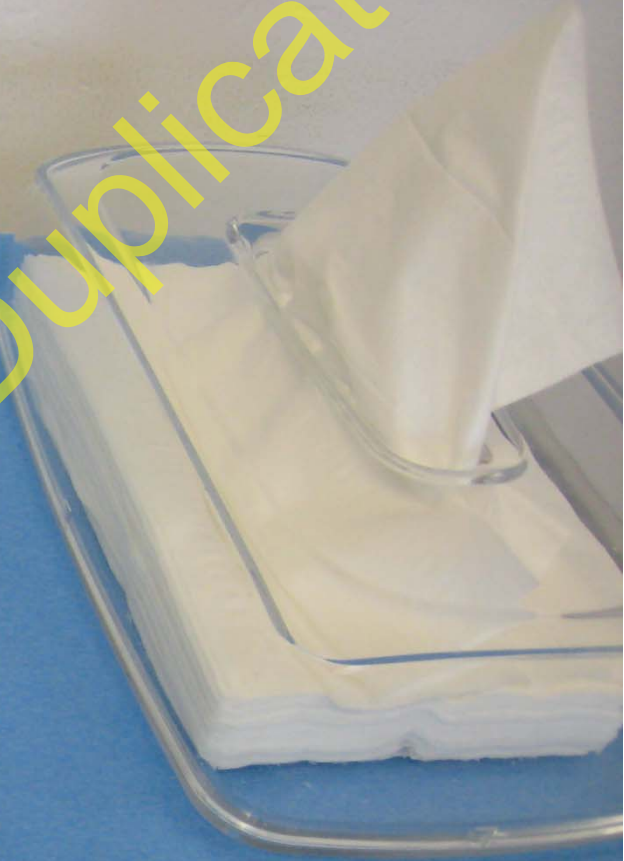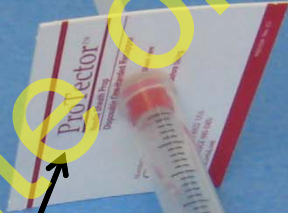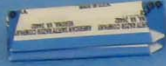

STICK SHIELD



Do Not Distribute or Duplicate

Stick Shield



Do Not Distribute or Duplicate

Sterilization Options:

Steam



Do Not Distribute or Duplicate

Sterilization Options:

Steam

Ethylene Oxide



Do Not Distribute or Duplicate




Do Not Distribute or Duplicate

Bishop-Harmon

**Delicate and Expensive Instrument:
Wrap and Sterilize Separately**



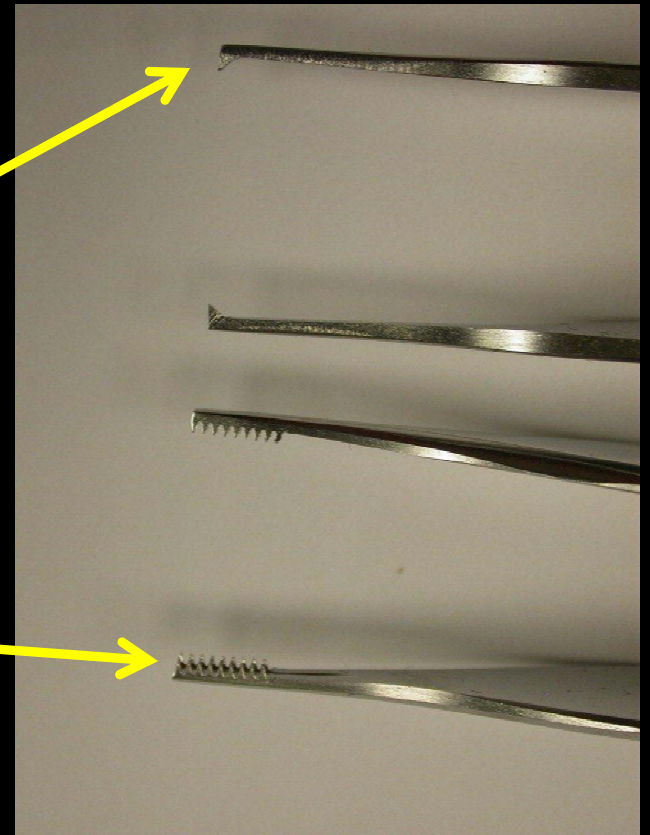

Brown Adson 1x2 vs. Multi-Tooth



This Stapler Can Close with Tension



**I HAVE NO FINANCIAL
INTEREST IN THIS COMPANY**



This Stapler Can Close with Tension



Backhaus

I HAVE NO FINANCIAL
INTEREST IN THIS COMPANY

Do Not Distribute or Duplicate

Periosteal Elevator

Nail Removal



Do Not Distribute Or Duplicate

Periosteal Elevator

Nail Removal

Blunt Excision (KA and Warts)



Periosteal Elevator



Nail Removal


Blunt Excision (KA and Warts)

Remove Periosteum

Do Not Distribute Or Duplicate


Webster Needle Driver with Smooth Jaws

40 Suture or Smaller

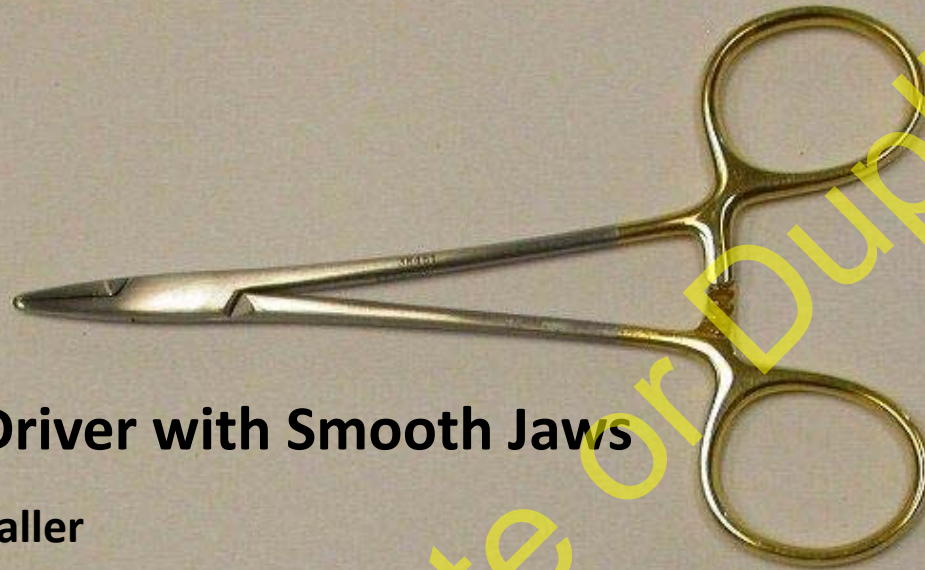

Crilewood with Serrated Jaws

30 Suture or Larger

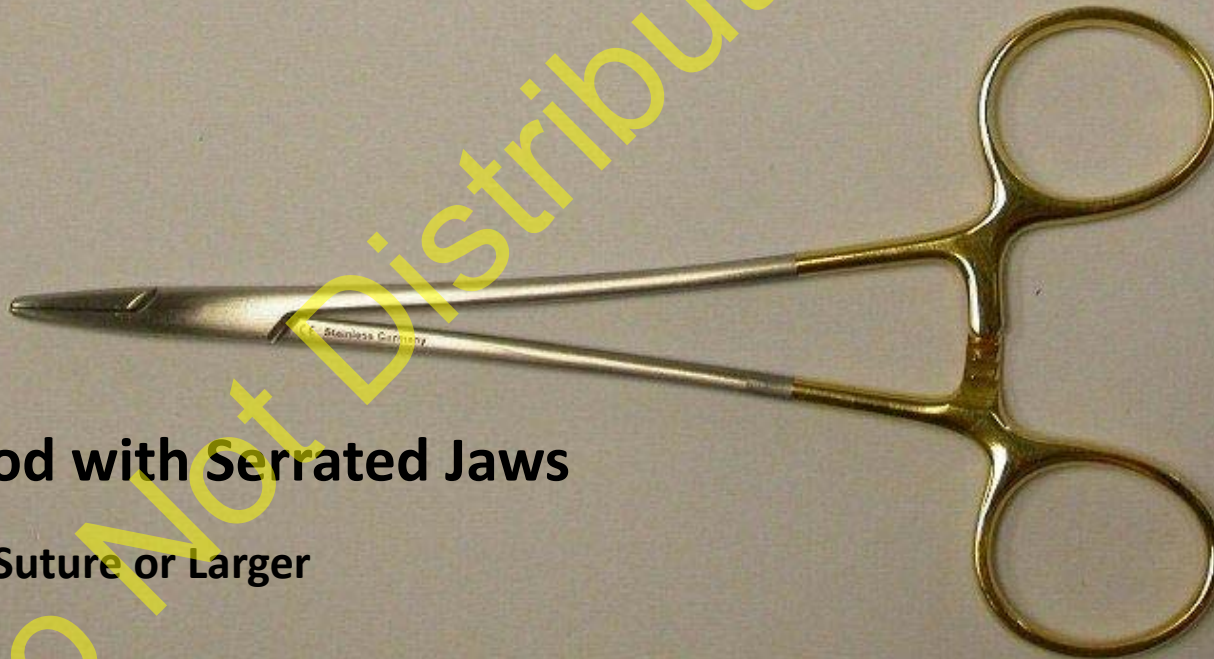

Locking Castroviejo Suturing Driver


Locking or Non-locking

5 $\bar{0}$ Suture or Smaller

Non-Locking Harder to Use (Mainly for Ocular Plastics)



Locking Castroviejo Suturing Driver

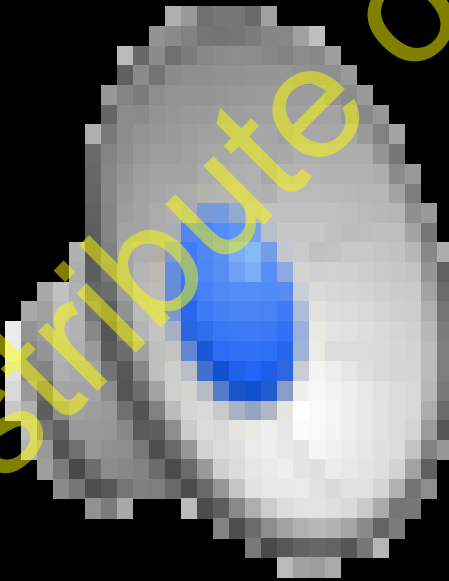

Do Not Distribute or Duplicate

Bausch + Lomb Magna Visor is more comfortable, and more durable.



Bausch + Lomb Magna Visor

Optivisor

**I HAVE NO FINANCIAL INTEREST
IN EITHER COMPANY**

Jeweler Forceps



Splinter Forceps



Do Not Distribute Or Duplicate

Tip is too Sharp and It Sticks to Tissue

Jeweler Forceps



Splinter Forceps



Do Not Distribute or Duplicate


Halstead Mosquito



Halstead Mosquito




Straight Version

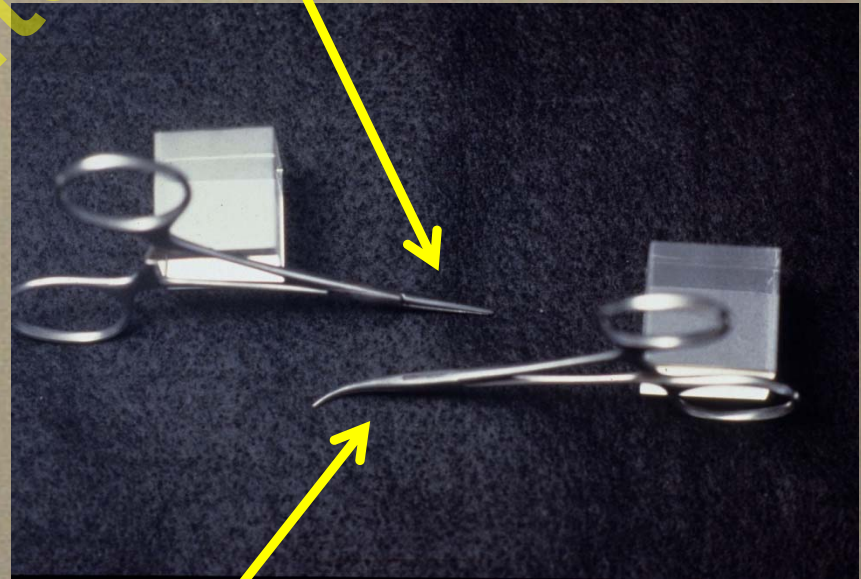

Curved Version

Do Not Distribute or Duplicate


Siegel Handle with Ruler

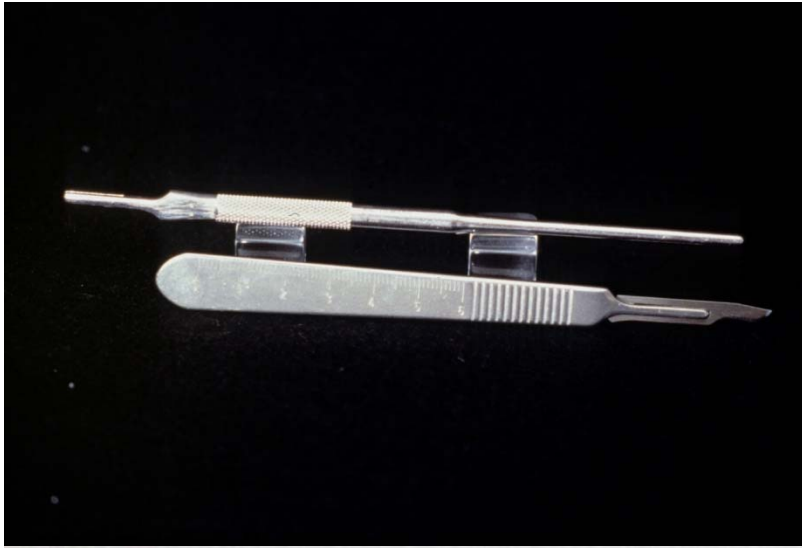

Siegel Handles are love or hate:
try **ONE** and decide.

Do Not Distribute or Duplicate

Siegel Handle with Ruler

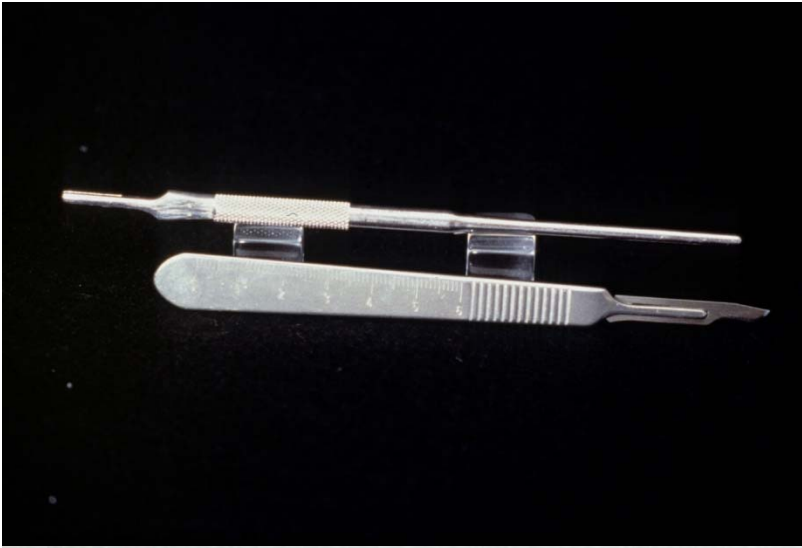

Siegel Handles are love or hate:
try one and decide.

I don't like them, but
Howard Steinman loves them.

All Knife Handles Should Have Rulers




Do Not Distribute or Duplicate

Re-Sterilizable Ruler In All Packs

Throw-Away Pads of Paper Rulers
Now Available




Roth Hook: No Better Than Any Other

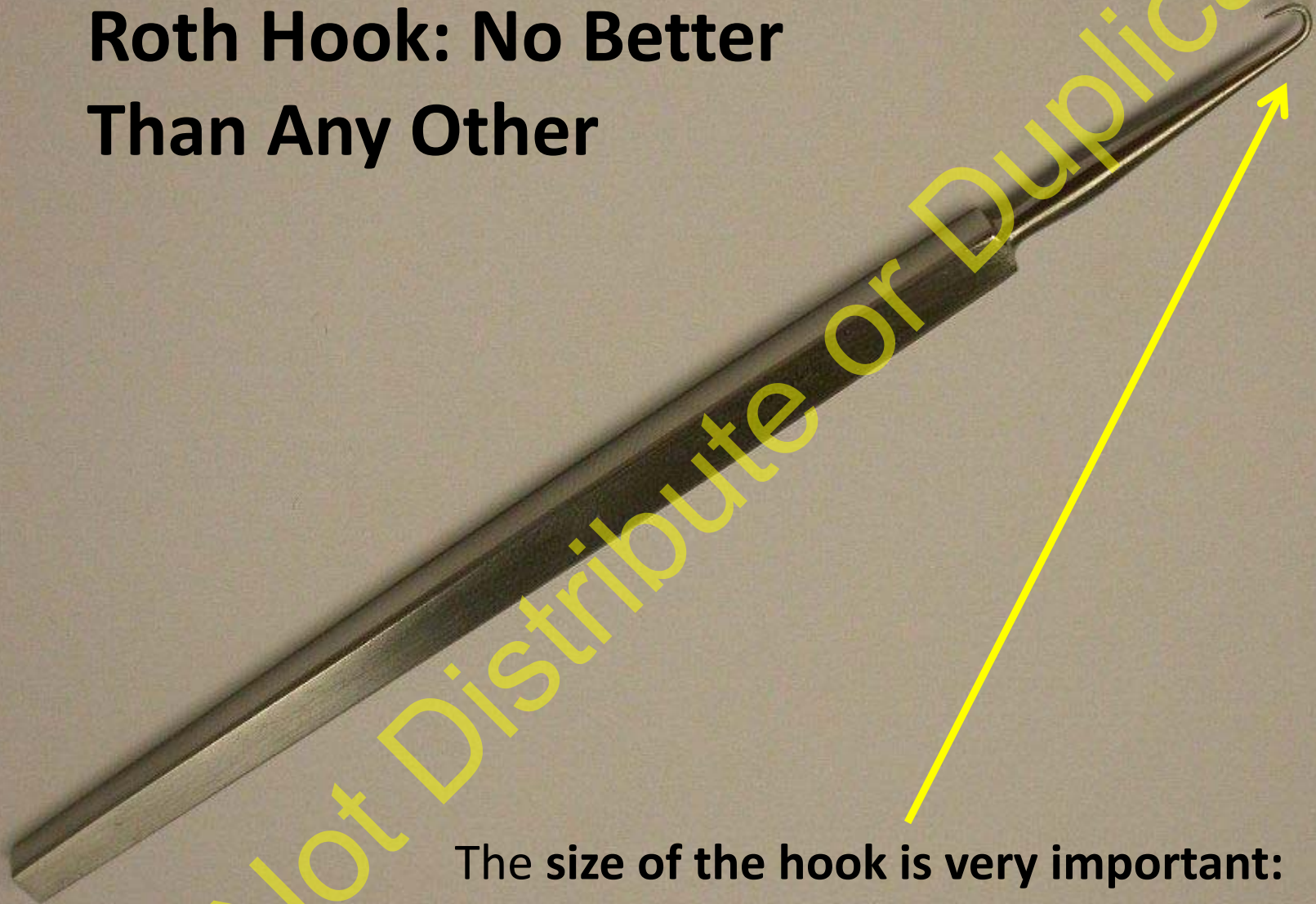

The size of the hook is very important:

Roth Hook: No Better Than Any Other

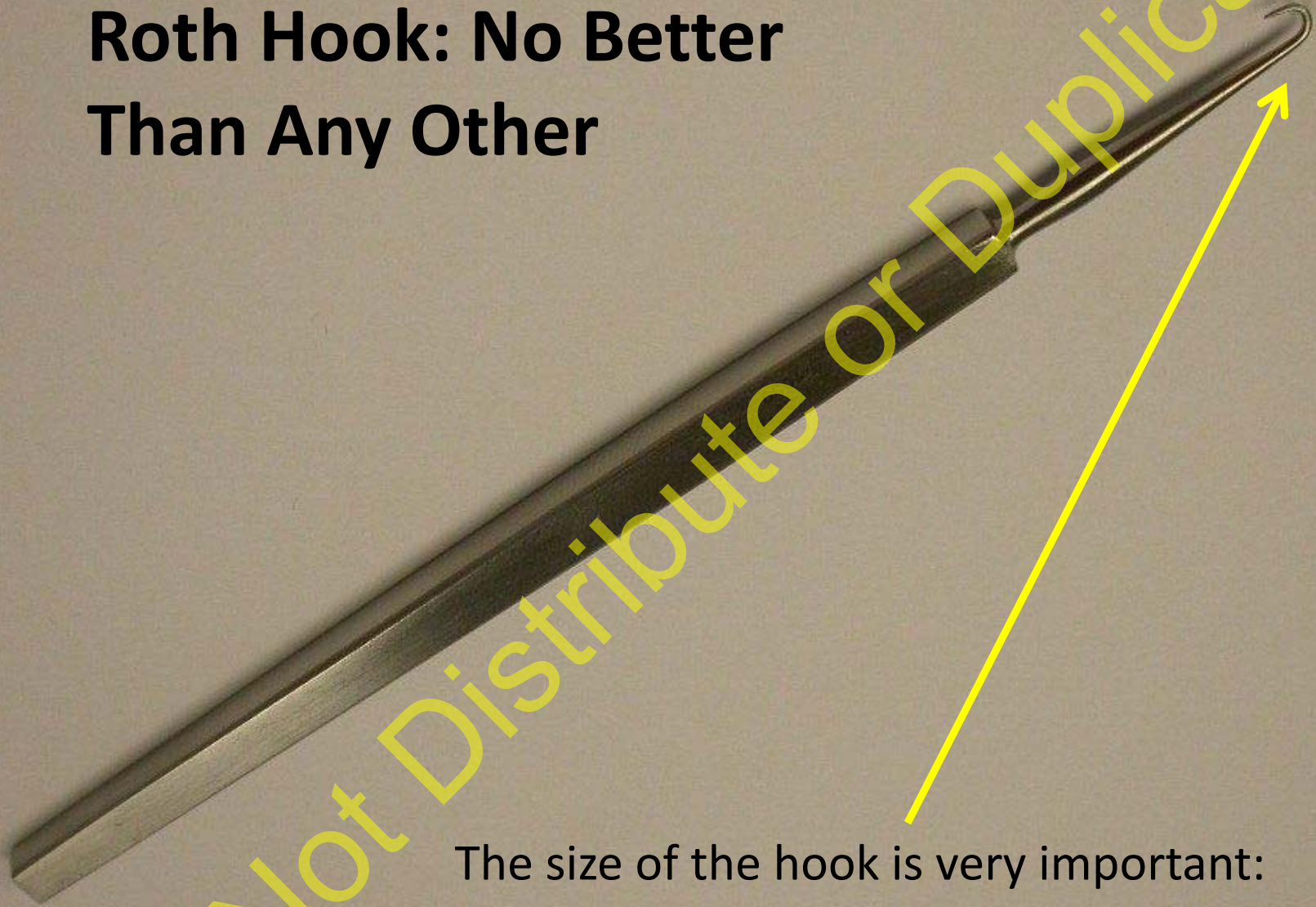

The size of the hook is very important:
Too **small** is **ineffective**,
Too **large** is **cumbersome**

Beaver Blade and Knife Handle System



Expensive and Seldom Used,
BUT Essential for Certain Cases

**Shown Here is the Blade I Use to
use for Tumors Extending into
the Ear Canal –available from
Tiemann**



Iris Scissors Are My Work Horse Undermining Scissors

Angled Scissors for
Medial Canthus and
Conchal Bowl





Iris Scissors Are My Work Horse Undermining Scissors

Angled Scissors for
Medial Canthus and
Conchal Bowl

Iris Scissors




I Use My Hammer and Chisels Surprisingly Often

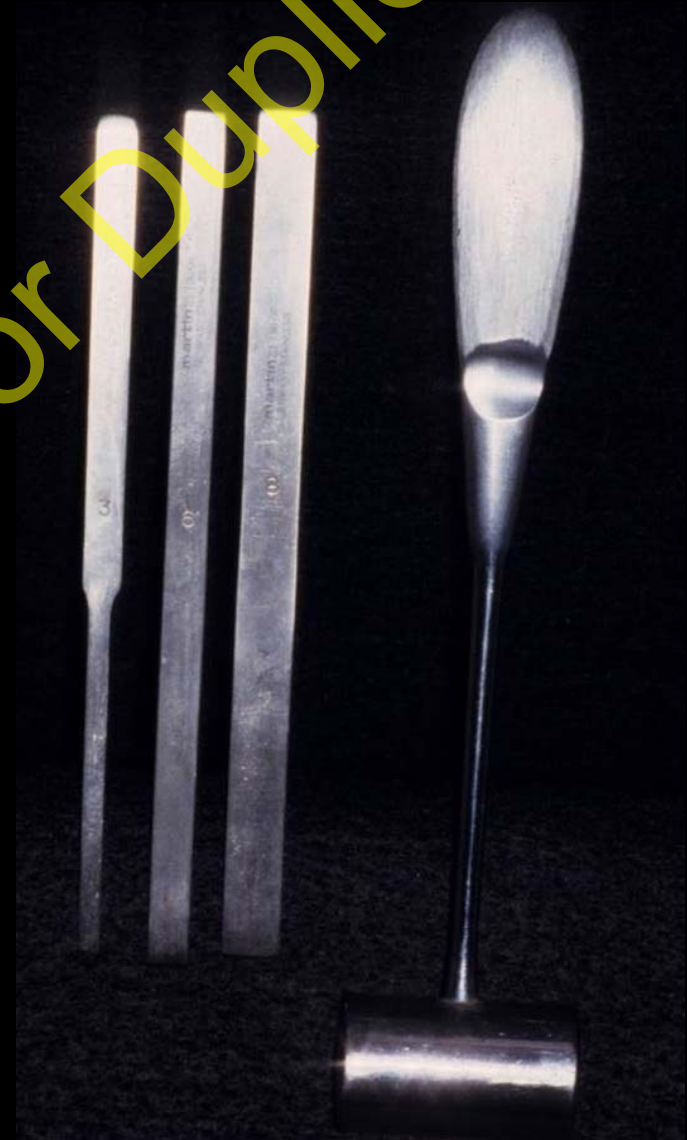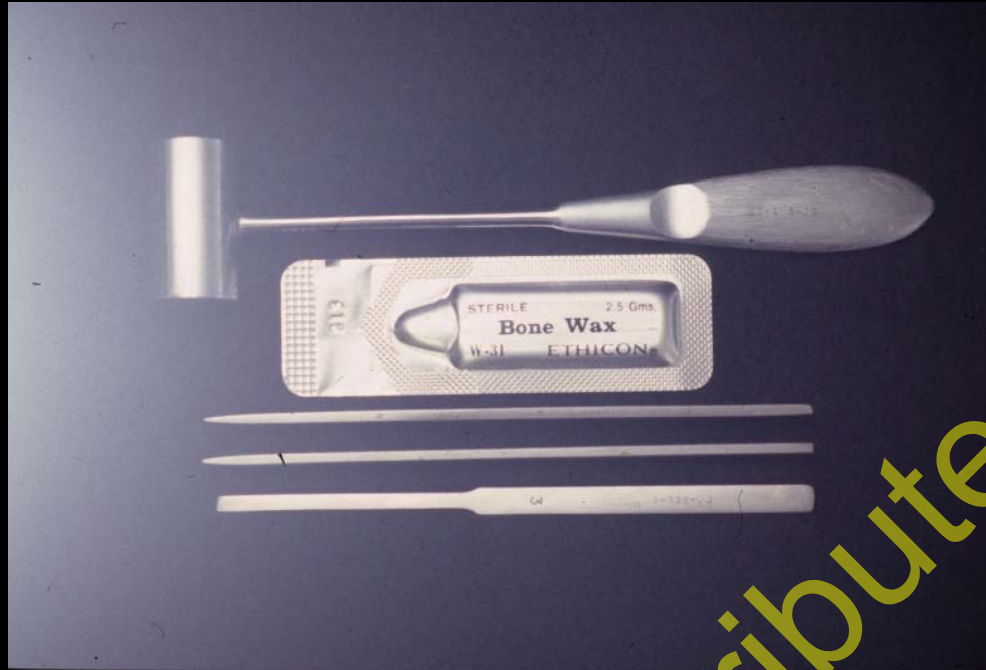


Do Not Distribute Or Duplicate

I Use My Hammer and Chisels Surprisingly Often

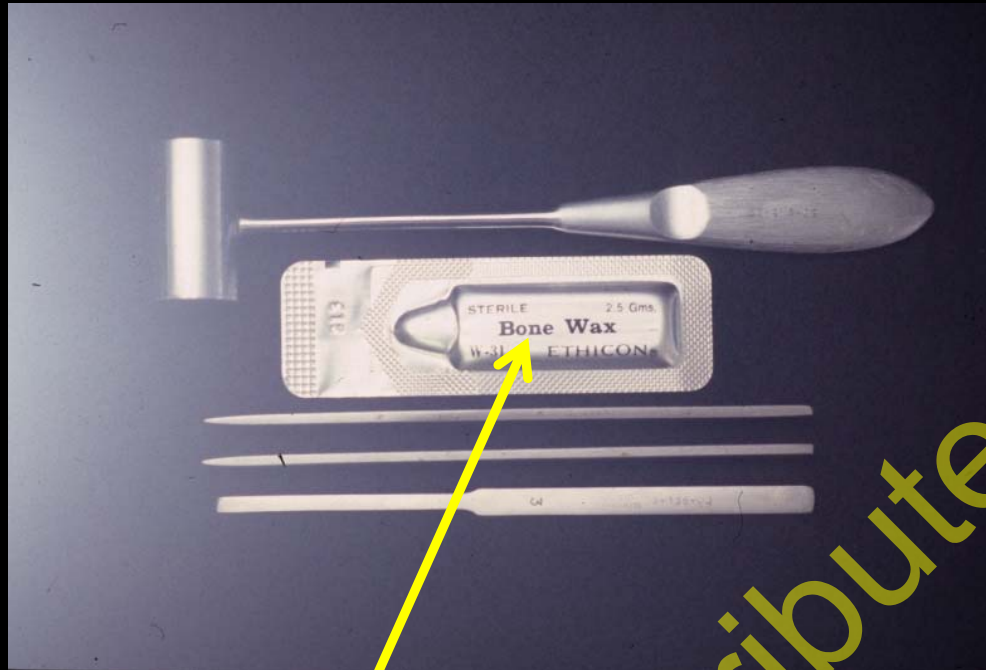

Don't Leave Home
Without Bone Wax

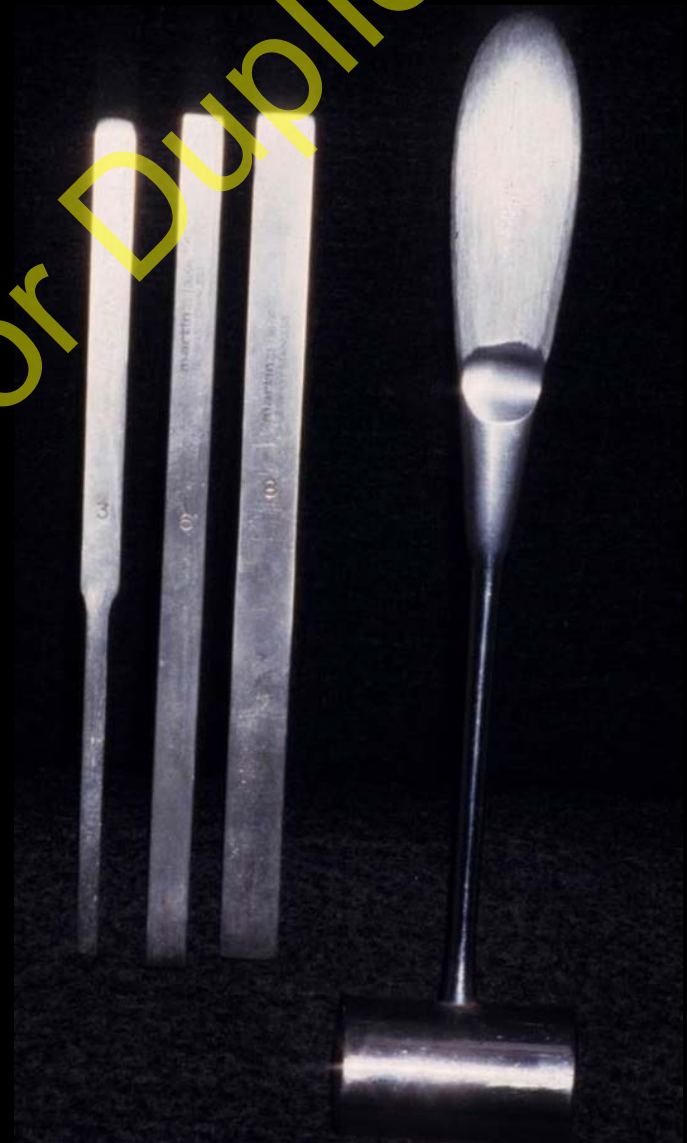

Do Not Distribute Or Duplicate

Chalazion Clamps



I Can't Seem to Make These Work Effectively
On Eyelids but May Be of Use in Other Areas

NAIL PACK


Spatula Nail Grabber



NAIL PACK

Spatula Nail Grabber

← Platypus Nail Grabber

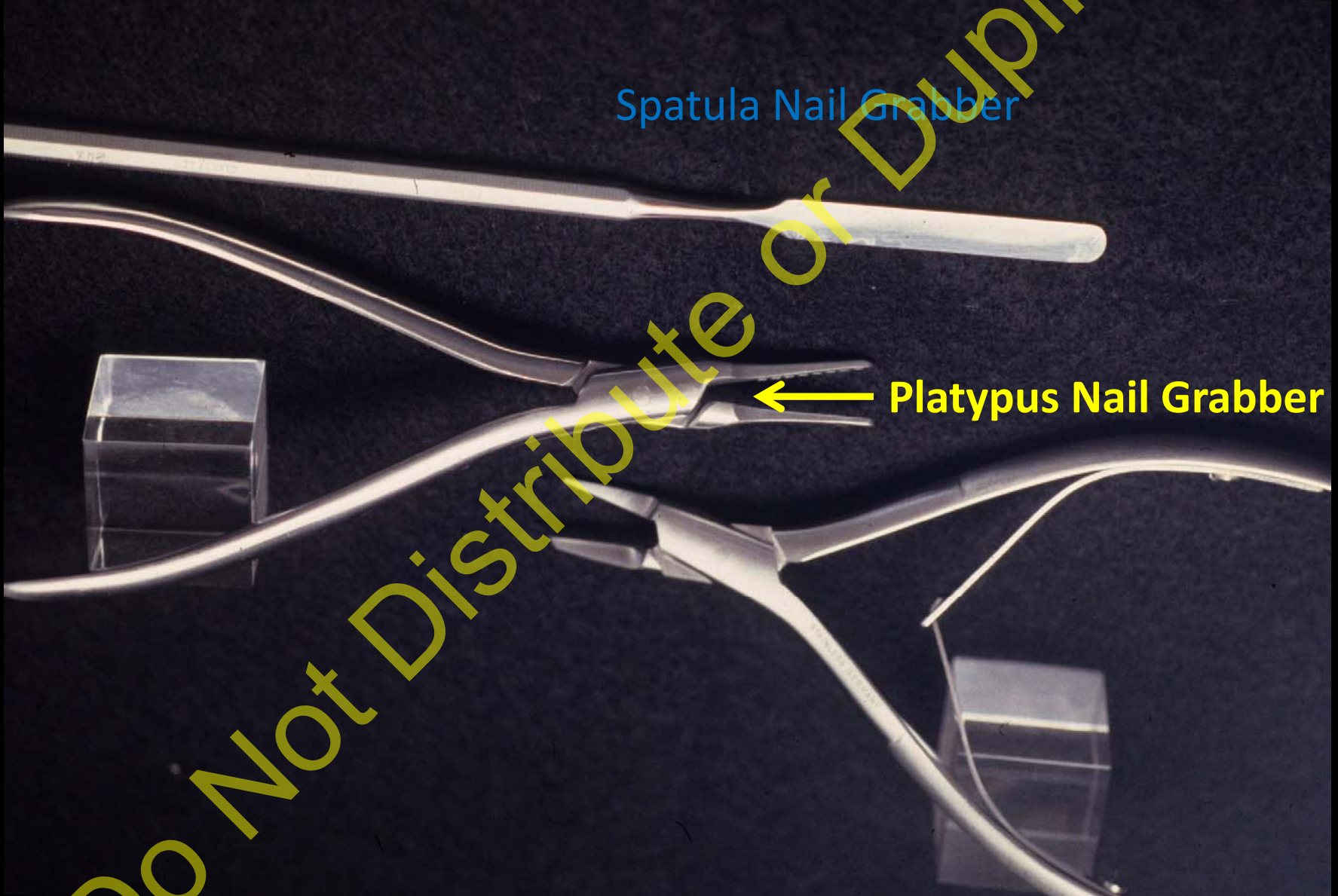


Do Not Distribute or Duplicate

NAIL PACK

Spatula Nail Grabber

Platypus Nail Grabber

Anvil Nail Splitter

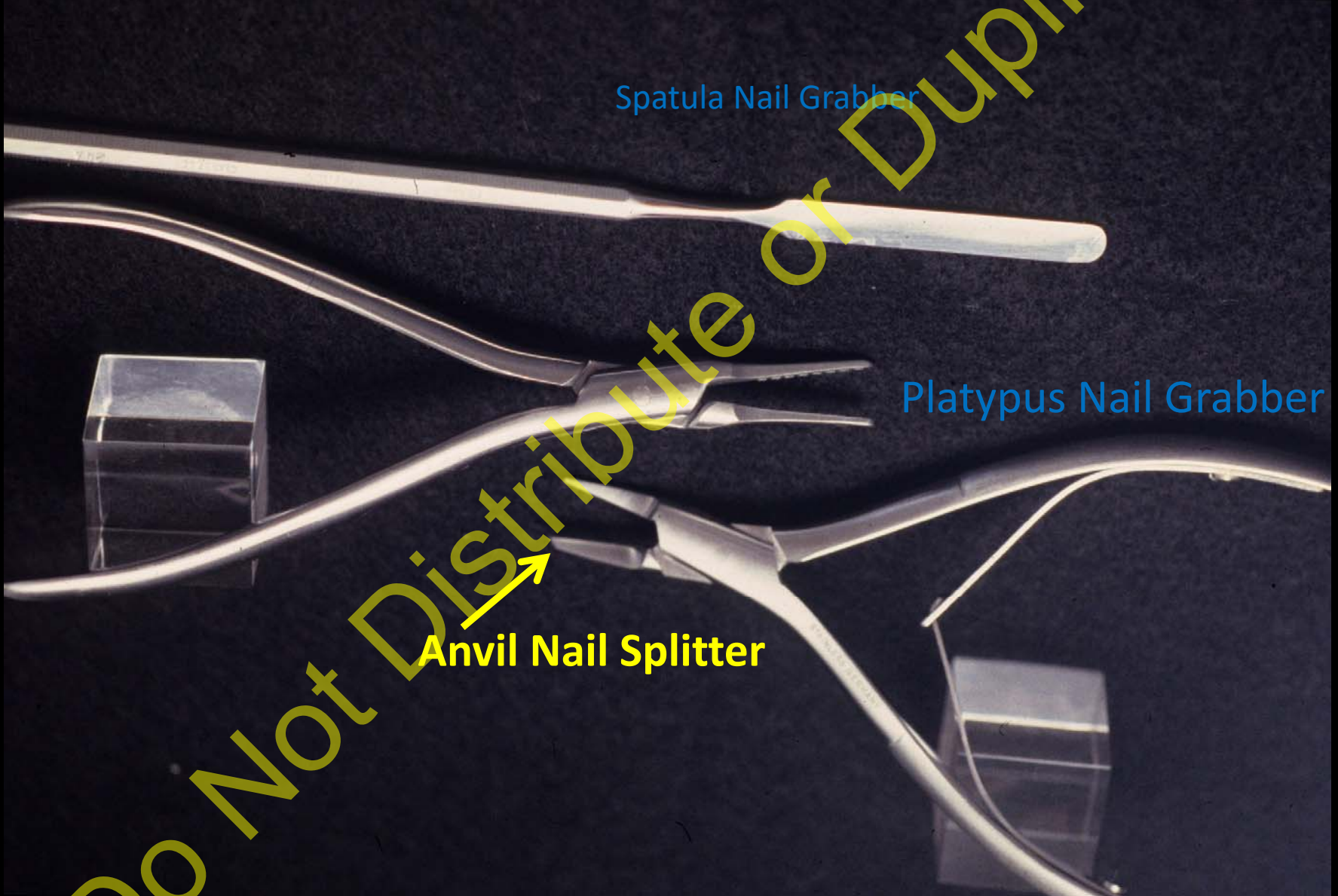


NAIL PACK


Spatula Nail Grabber

Platypus Nail Grabber

Anvil Nail Splitter


Periosteal Elevator

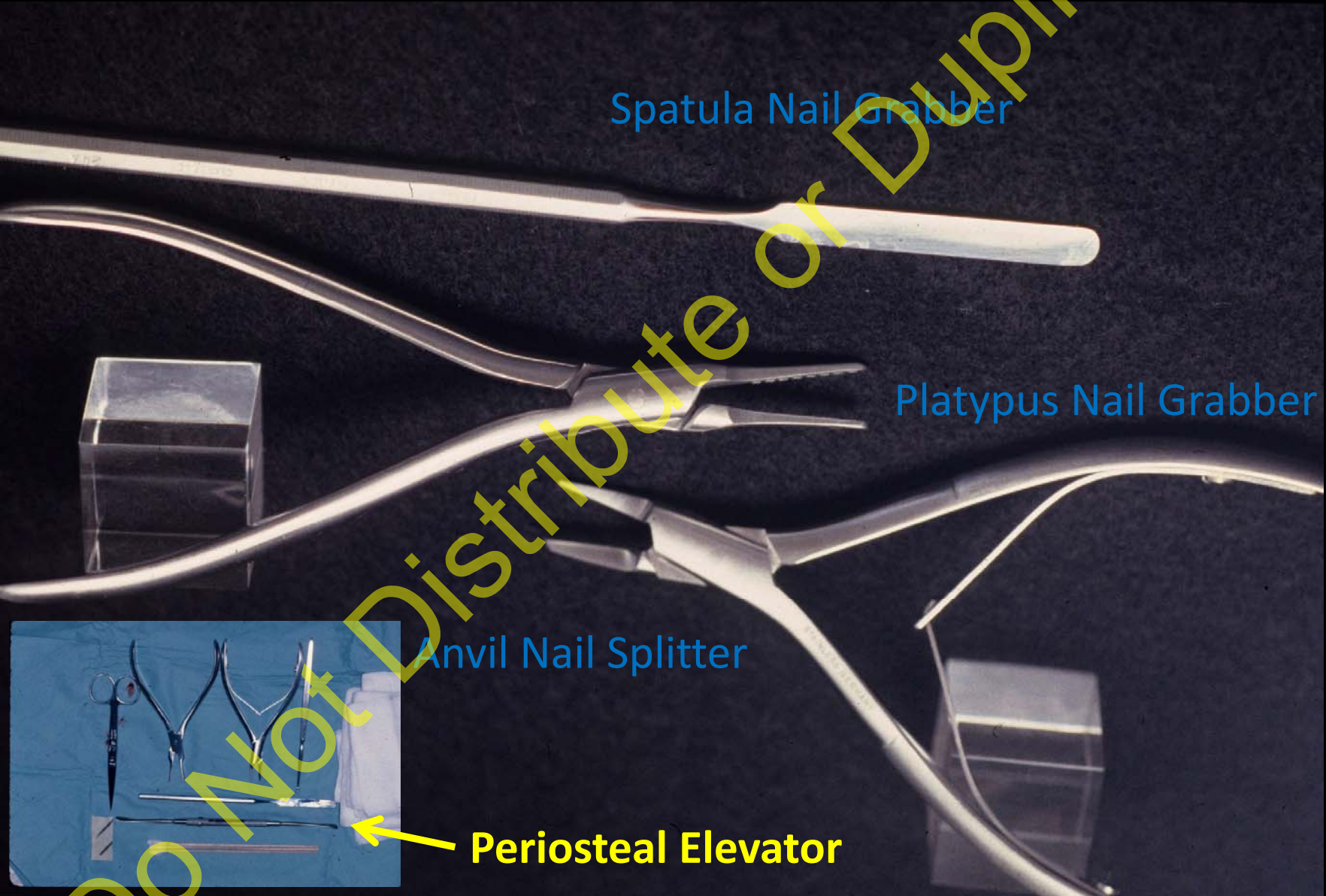

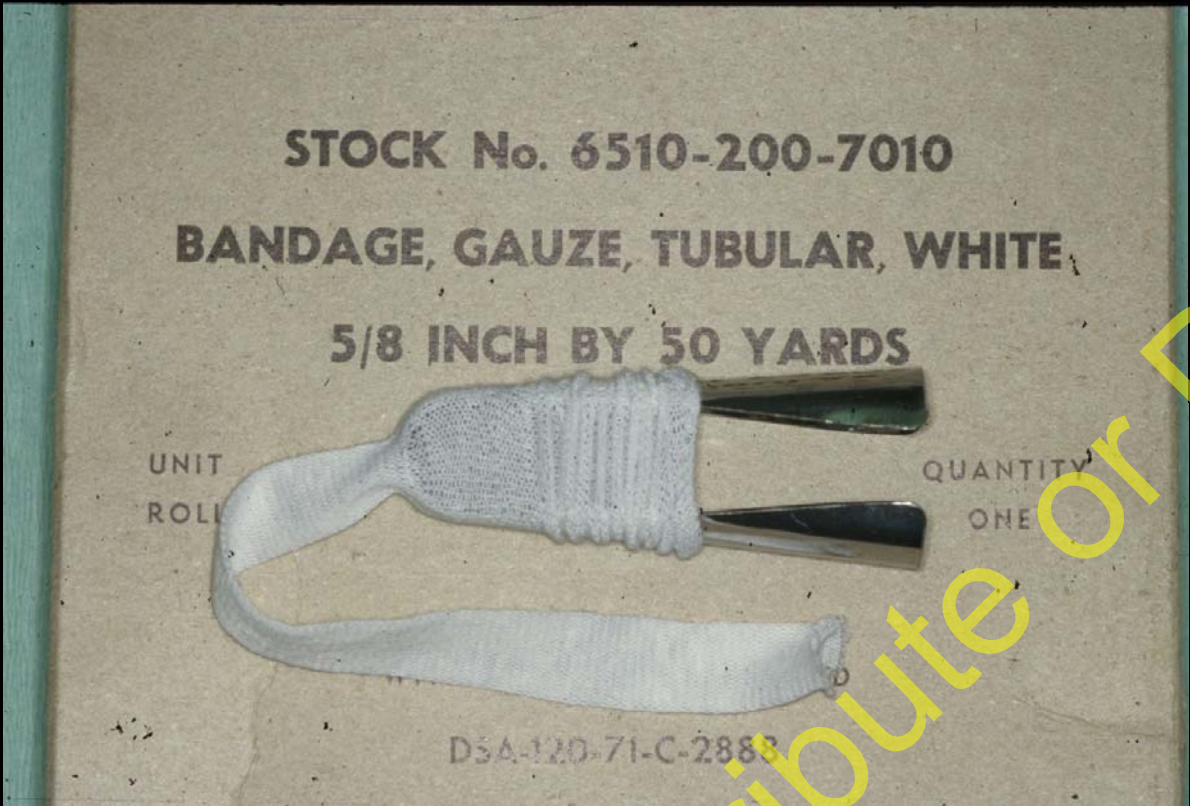

Tube Gauze Comes Small
Enough Sizes For A Single Finger



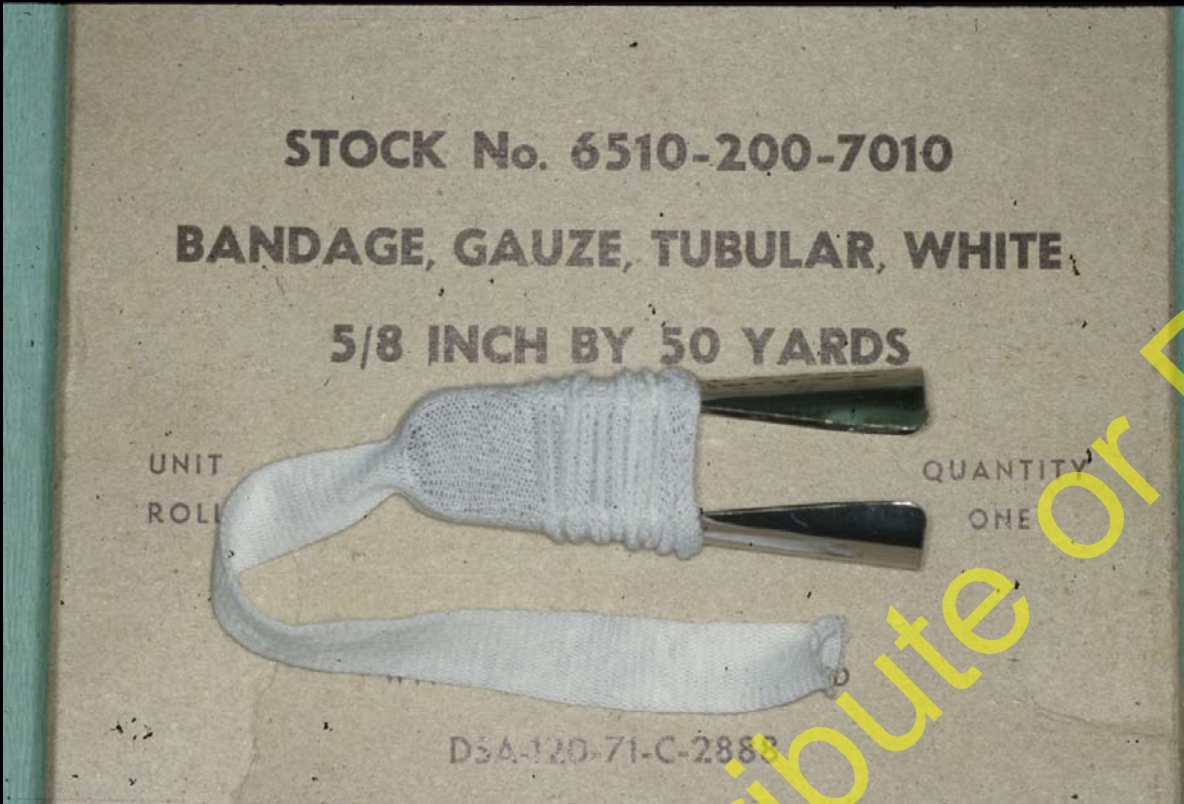

Tube Gauze Comes Small Enough Sizes For A Single Finger



Tube Gauze Also Comes Big Enough Sizes For An Entire Head

Do Not Distribute or Duplicate

Gradle Ocular Scissors



Alligator Hair Clamp



Do Not Distribute or Duplicate

Sterile Sleeve



Towel Clamp Holds Bovie in Place